addition mini Junior/Sensitive

GB
Instructions for use

Precision impression series for babies, infants & adolescents

Specially designed for the needs of children and the specific requirements of paediatric adaptions.

Junior Children and adolescents, easy-flow, neon-green with scent of gummi bear

Sensitive Babies, infants and sensitive patients, ultra-easy flow, does not contain dyes, pigments, flavouring or preserving agents

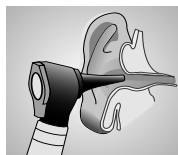
Fig. 1

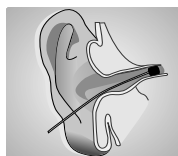
Fig. 2

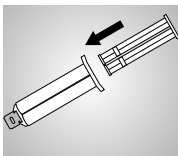
Fig. 3

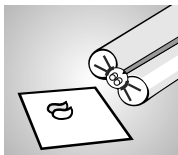
Fig. 4

Fig. 5

Fig. 6

Fig. 7

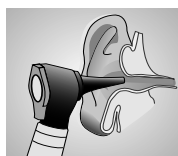
Fig. 8

1. Ear preparation

Before taking an impression, the auditory canal and the eardrum have to be examined (Fig. 1). If distinctive features are found (e.g. inflammations or perforation of the eardrum, tubes) impression must not be taken. The auditory canal has to be cleaned. In case of insufficient cleaning, remaining cerumen can cause retarded setting resp. inhibition of setting process.

Before taking the impression, the eardrum must always be protected with an impression plug/cotton placed at the end of the external auditory canal in front of the eardrum (Fig. 2). The use of **vented impression pads** is recommended. Now **addition mini Junior/Sensitive** is injected directly into the prepared ear of the patient as described under point 2.

2. Mixing and applying

Dispensing is performed with the minimix system. Insert the plunger in the body of the cartridge (Abb. 3). If necessary, the mini mixing gun can be used. Remove the cap from the cartridge by turning. Carefully extrude a small amount of material to ensure proper flow from both orifices (Fig. 4). Insert mixing cannula into guide grooves on cartridge and turn clockwise to lock in position (Fig. 5). If required, a piece of silicone tubing (Flexi tube) can also be pushed onto the mixing cannula. Again, bleed small amount to ensure uniform and even material flow. Leave mixing cannula on cartridge after use, thus sealing the cartridge. Press out a small amount and discard.

In order to allow a true rendition of the ear canal without deformation, the impression material must have a low viscosity and be pourable. The higher the viscosity, the higher the pressure build-up and thus the deformation of the child's ear canal. The good flow properties of **addition mini Junior/Sensitive** are therefore desirable and cannot be compared with those of conventional materials for adults.

3. Dispensing

addition mini Junior/Sensitive flows ultra-easily, thus allowing pressureless, constant filling. Slowly apply the impression material into the auditory canal until the concha is filled completely (Fig. 6). Keep the end of the cannula/tube immersed in the material during dispensing. Tip: Hold the head slightly at an angle or put on some foil in order to prevent any oozing out. Leave mixing cannula on cartridge after use, thus sealing the cartridge. After curing remove the impression carefully from the ear (Fig. 7). Finally it is essential to re-examine the ear to make sure it has not been damaged (Fig. 8). The level of the minimix cartridge can be checked according to how far the plunger is protruding.

Important working hints

- Cured impression materials are chemically inert, pots on clothing should be avoided.
- Be sure to remove all traces of impression material from the patient's ear.
- Latex gloves and latex contaminated surfaces as well as cerumen, cremes and resins may inhibit the setting reaction of **addition mini Junior/Sensitive** (we recommend standard gloves made of nitrile rubber or polyethylene).
- A careful preparation (point 1) as well as a thorough cleaning of the auditory canal is required before taking an impression.
- For processing soft, high-elasticity silicone materials with a low final hardness (< 30 Shore A), special grinding sleeves are generally better than conventional metal milling machines. We recommend **Detax grinding sleeves**: Round, Ø 5 mm (Art.No. 03006, 10 pcs.) or Pointed, Ø 5 mm (Art.No. 03011, 10 pcs.).

Safety advice

The impression taking of the external auditory canal has to be executed by trained specialist staff only. The working instructions and precautionary measures have to be strictly observed. DETAX is not liable for any damage caused by improper application of the impression material.

Further information: Silicone based materials are proven a million times. On condition of a proper application, undesired effects are not to be expected. However, reactions of the immune system like allergies, irritations, cannot be absolutely excluded. In case of doubt, we recommend to make an allergy test before the application of the material.

Indications for use:

- Ear impression with:
 - Babies, infants, adolescents
 - Sensitive patients

Technical data:

- **Mixing volume:** 10 ml (minimix cartridge)
- **Mixing ratio:** 1:1
- **Colour code:**
 - Junior** Base: transparent catalyst: neon-green
 - Sensitive** Base: transparent catalyst: white
- **Application time:** approx. 30 sec.
- **Time in the ear:** approx. 2:00 min. at body temperature
- **Total setting time:** approx. 2:30 min.
- **Final hardness:** approx. 25 Shore A
- **Strain in compression:** approx. 10 % ± 1 %
- **Recovery from deformation:** 99,8 %
- **Linear dimensional change:** < 0,2 %
- **Application:** At 23 °C ± 2 °C / 73 °F ± 4 °F, 50 ± 5 % rel. humidity
- **Storage:**

25°C

77°F

15°C

59°F

Ordering information:

addition mini Junior
Standard packing **03439**
4 x 10 ml minimix cartridges
8 mixing cannulas
minimix, green
1 plunger, flowchart

Refill bag 03440
8 x 10 ml minimix cartridges
1 plunger, flowchart

addition mini Sensitive
Standard packing **03436**
4 x 10 ml minimix cartridges
8 mixing cannulas
minimix, green
1 plunger, flowchart

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